

Central Lancashire Online Knowledge (CLoK)

Title	Lipid droplet detection by the cavity perturbation method
Type	Article
URL	https://clok.uclan.ac.uk/id/eprint/5381/
DOI	https://doi.org/10.1088/1742-6596/307/1/012021
Date	2011
Citation	Blakey, R T, Mason, A, Rolph, C E, Bond, Gary and Al-Shamma'a, A I (2011) Lipid droplet detection by the cavity perturbation method. Journal of Physics: Conference Series, 307 (1). 012021. ISSN 1742-6596
Creators	Blakey, R T, Mason, A, Rolph, C E, Bond, Gary and Al-Shamma'a, A I

It is advisable to refer to the publisher's version if you intend to cite from the work. https://doi.org/10.1088/1742-6596/307/1/012021

For information about Research at UCLan please go to http://www.uclan.ac.uk/research/

All outputs in CLoK are protected by Intellectual Property Rights law, including Copyright law. Copyright, IPR and Moral Rights for the works on this site are retained by the individual authors and/or other copyright owners. Terms and conditions for use of this material are defined in the http://clok.uclan.ac.uk/policies/



Home Search Collections Journals About Contact us My IOPscience

Lipid droplet detection by the cavity perturbation method

This article has been downloaded from IOPscience. Please scroll down to see the full text article.

2011 J. Phys.: Conf. Ser. 307 012021

(http://iopscience.iop.org/1742-6596/307/1/012021)

View the table of contents for this issue, or go to the journal homepage for more

Download details:

IP Address: 193.61.255.75

The article was downloaded on 02/07/2012 at 16:18

Please note that terms and conditions apply.

doi:10.1088/1742-6596/307/1/012021

Lipid droplet detection by the cavity perturbation method

R T Blakey¹, A Mason¹, C E Rolph², G Bond³ and A I Al-Shamma'a¹

Email: r.t.blakey@2010.ljmu.ac.uk

Abstract. There are currently no point-of-care diagnosis strategies available to indicate the presence of neoplasmic growth. This research aims to develop a novel diagnostic strategy based on detecting TAG accumulation in cells. This element of the research is a preliminary experiment to prove the concept of detecting TAG lipid droplets in YEPD media. It was found that a change in mono-unsaturated concentration can be detected by the frequency shift in a resonant cavity. The dielectric constant of TAG vegetable oils was calculated at 2.34-2.39. It was also found that concentrations of lipid droplet can be differentiated up to 5% (v/v).

1. Introduction

Despite decades of research, cancers or malignant neoplasms continue to significantly contribute to patient morbidity and mortality. Neoplasms occur when cell proliferation exceeds cell apoptosis (naturally initiated cell death) which results in the formation of a tumour or increased blood cell volume depending on the origin of the cancerous cell. Tumours can be problematic to detect in the early stages due to symptoms not arising until later stages of development. When tumours are detected in earlier stages treatments become less invasive, simpler and cheaper resulting in increased patient comfort and prognosis. Screening programs have been extremely successful in lowering the incidence of late stage cancers (e.g. breast [1], cervical [2]) but are limited to a specific type and site of cancer also requiring specialised facilities, staff and transport (i.e. laboratories). Currently there are no quantitative or qualitative point-of-care diagnostic assays available to indicate neoplasmic tissue. This communication proposes a novel diagnostic strategy that can detect a tumour in the body from any origin.

This is proposed through the measurement of triacylglycerol (TAG) stores within cancerous cells and specialised cells of the immune system called dendritic cells. Under normal circumstances, most cells store small amounts of TAGs in organelles called lipid droplets (lipid bodies, adiposomes) and are used to synthesise lipid membranes, other molecules of the lipid family and as a source of

¹ School of Built Environment, Liverpool John Moores University, Liverpool L3 3AF, UK

² School of Pharmacy and Biomedical Sciences, University of Central Lancashire, Preston PR1 2HE, UK

³ School of Forensic and Investigative Sciences, University of Central Lancashire, Preston PR1 2HE, UK

r.t.blakey@2010.ljmu.ac.uk

doi:10.1088/1742-6596/307/1/012021

chemical energy. However, due to cellular mutations and/or abnormal cell to cell signalling, the metabolism of the cell is altered resulting in TAGs being accumulated in such lipid droplets. Depending on the mutations that lead to the neoplasm, a lipid droplet can constitute a significant sum of cell volume [3]. Hence, the cell concentration of TAG is an excellent indicator or biomarker of cell metabolism and proliferative state. Furthermore, excess TAGs released into circulation are "mopped up" by cells of the immune system. Dendritic cells are specialised cells of the immune system which scavenge for non-self or otherwise foreign molecules and present digested material to other cells of the immune system (T-cells). Recent research has found that dendritic cells from tumour bearers accumulate TAG due to exacerbated macrophage scavenger receptor 1(MSR1) mediated uptake [4]. This relationship is complex and it is unclear whether TAG accumulation inhibits the immunogenic actions of dendritic cells facilitating the growth of tumours or the result of a tumour being present. Dendritic cells have the advantage of being present in the blood stream while migrating to specific sites around the body allowing cells to be easily isolated from the body (e.g. blood sample) to be characterised, removing the need for invasive diagnostic strategies such as surgery or biopsy. Lipid droplet accumulation also occurs in a range of cells specifically related to conditions of infection and inflammation [3].

Many different impedance or dielectric spectroscopy techniques have been developed in an array of cancer diagnostic applications and show promising results. Approaches vary depending on the site and type of cancer from characterising isolated individual cells to non-invasive *in-vivo* imaging techniques [5-7]. Dielectric spectroscopy analyses of cellular material are based on dielectric properties at different frequencies. Biological cells exhibit a characteristic decrease in permittivity as a function of increasing frequency known as α , β , δ and γ dispersions [8]. The aforementioned review is particularly useful to the biologist and engineer alike who may not be familiar with electromagnetic biological cell interactions.

The following experiment provides a simple preliminary proof of concept trial to determine if different conformations of TAGs and different concentrations of TAG droplets can be differentiated based on their dielectric properties at microwave frequencies. This is calculated using the small perturbation method whereby a sample introduced to a resonating cavity shifts the frequency of resonant modes and decreases the quality factor dependant on the dielectric characteristics of the sample under test [9-11]. Resonant based sensors have proved to be extremely efficient in detecting the small changes in the biochemical make up and biological characteristics of a cell which can be calculated based on the measured electromagnetic properties of the cell sample [12].

2. Methods

2.1. Chemicals and reagents

Yeast extract, bacto-peptone, dextrose (YEPD) were purchased from Sigma-Aldrich UK. Grapeseed, Groundnut, Olive and Sunflower Oils were purchased as non-branded products from a well known U.K. supermarket.

2.2. Preparation of TAG and TAG lipid droplet samples

Samples of 100% grapeseed, groundnut, olive and sunflower oil were prepared in sample tubes ready for dielectric measurement. Secondly, olive oil was suspended in YEPD media to yield concentrations of 50, 10, 5, 1, and 0.5 (% v/v) oil. YEPD was prepared to the following protocol: 1 % (w/v) of yeast extract, 2 % (w/v) of bacto-peptone and 2 % (w/v) of dextrose diluted in distilled water. Before measurement, the tube to be sampled would be vortex until a homologous suspension was achieved.

2.3. Dielectric measurement and calculation

Measurements were performed using an Agilent Technologies (Hewlett Packard) 8720 ET Vector Network Analyser (VNA). The instrument was set to generate a signal between 1.5-1.9 GHz over 1601 data points and calculate the S_{21} (transmission) parameters. Measurements were made through a

doi:10.1088/1742-6596/307/1/012021

custom cylindrical resonating cavity, designed and constructed in house. The cavity was designed to resonate at 1.7 GHz at mode TM_{010} .

The dielectric constant, dielectric loss and loss tangent of the sample were then calculated using the following equations [11].

$$\varepsilon' = \frac{V_c (f_c - f_s)}{2V_s f_s} + 1$$

$$\varepsilon'' = \frac{V_c}{4V_s} \left(\frac{1}{Q_s} - \frac{1}{Q_c}\right)$$

$$\tan \delta = \frac{\varepsilon''}{\varepsilon'}$$
(1)

(3)

where ε' and ε'' is the dielectric constant and dielectric loss respectively, f_c and f_s are the resonant frequencies of the empty cavity and loaded cavity respectively, V_s and V_c are the volumes of the sample and cavity respectively and Q_c and Q_s are the quality factors of the empty cavity and the loaded cavity respectively.

3. Results

Figure 1 shows the S₂₁ parameters for the vegetable oil standards and empty cavity as a function of frequency illustrating frequency perturbation for groundnut oil of 17.1 MHz, olive oil of 17.4 MHz, sunflower oil of 17.5 MHz and grapeseed oil by 17.7 MHz to 1.69 GHz.

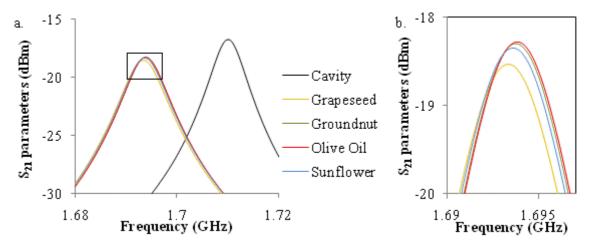


Figure 1a. S₂₁ parameters of TAG samples (magnified in figure 1b.) and the unloaded cavity. S₂₁ parameter frequency spans from 1.68 to 1.72 GHz.

Figure 2 shows the S₂₁ parameters of differing concentrations of Olive Oil suspended in YEPD. It can be seen that as the concentration of olive oil decreases, the resonant frequency of the sample decreases or shifts to the left away from the resonant frequency of the empty cavity. Results can be seen in table 1. The total frequency shift between the empty cavity and 0.5% oil is 135 MHz to 1.576 GHz (Peak E). Peak D and E consist of two curves on the below graph, 10% -5% and 1%-0.5% respectively. The quality factor of the detected parameters decreases as YEPD media is introduced into the sample.

doi:10.1088/1742-6596/307/1/012021

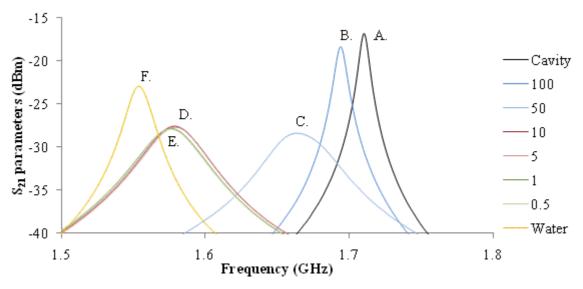


Figure 2. S_{21} parameters of lipid droplets present in differing concentrations (olive oil % v/v). Peaks from right to left are A: Cavity. B: 100%. C: 50%. D: 10% and 5%. E: 1% and 0.5%. F: Water.

3.1 Dielectric calculations

The dielectric constant, dielectric loss and loss tangent of the samples were calculated using (1-3). The same parameters were also calculated for the TAG-YEPD suspensions. The results can be seen below in table 1. It can be seen that grapeseed and sunflower oil have the greatest shift in frequency resulting in a higher dielectric constant. This is due to the levels of unsaturated acyl chains present on TAG molecules that make up the oil. Groundnut and Olive Oil are composed of approximately 53% (w/w) and 72% (w/w) monounsaturated fatty acids respectively as compared to 21% (w/w) and 16% (w/w) of sunflower and grapeseed oil. This changes the dielectric characteristics of the sample resulting in a frequency shift. This is in agreement with previous research that states the dielectric constant of TAGs increases as the level of unsaturated fatty acids (FA) increases [13]. The Q-factor of the oil samples decreases indicating an increase in dielectric loss.

Table 1. Calculated dielectric properties of oil samples and suspensions.

	Frequency shift (MHz)	Q-factor	Dielectric constant	Dielectric loss	Loss tangent
TAG					
Grapeseed	17.7	199.21	2.387	0.100	0.003
Groundnut	17.1	199.28	2.340	0.100	0.004
Olive	17.4	199.25	2.364	0.100	0.004
Sunflower	17.5	201.61	2.372	0.006	0.004
Olive Oil (% v/v)					
100	17.0	199.25	2.364	0.100	0.003
50	47.0	35.03	4.749	1.572	0.331
10	132.5	38.04	12.142	1.422	0.117
5	132.0	38.99	12.096	1.380	0.114
1	134.5	38.93	12.324	1.382	0.112
0.5	135.0	39.40	12.370	1.362	0.110

doi:10.1088/1742-6596/307/1/012021

The dielectric constant of the oil suspensions increases as the concentration of YEPD increases. This is due to the water content of the media shifting the frequency closer to that of the water standard. Dielectric losses also increase as YEPD is introduced. This is due to the ion content of the media resulting in energy being lost as dielectric heating.

4. Conclusions and Discussion

This experiment has proved that differing conformations of lipids and different concentrations of lipid droplets can be distinguished based on the electromagnetic properties of the suspensions which are apparent in the frequency shift of the samples at this frequency range. However, limitations include the sample size being too large for this method. Accuracy of the technique increases as the sample size decreases. The technique also has the drawback of being inaccurate on low loss samples [11]. The current resonator and sample volumes are sensitive enough to detect lipid droplet concentration up to 5% (v/v). Further modification of the resonator cavity used to detect the lipid droplets will positively advance the accuracy of the calculations involved with this method. The homogeneity of the sample may also have introduced errors. As the sample is introduced to the cavity the TAG droplets begin to converge into larger droplets until the TAG begins to separate into a different phase resting on the YEPD media. This results in the lipid droplet concentration being measured actually being lower than the newly vortex suspension. Stabilisation of the lipid droplets using phospholipids would improve the homogeneity of the suspensions improving measurement accuracy.

The next stage of the research will look to detect lipid droplets within living yeast cell cultures. A number of complex theoretical models exist that can calculate the size and shape of the cell and composition of the cytoplasm including organelles from the dielectric characteristics of a sample. Such models will prove to be crucial in the correct determination of TAG accumulation within living cells.

References

- [1] Wyld L and Ingham C E 2007 Screening of the population for breast cancer Surgery 25 254–6
- [2] Peto J, Gilham C, Fletcher O and Matthews F E 2004 The cervical cancer epidemic that screening has prevented in the UK *Lancet*. **364** 249–56
- [3] Bozza P T and Viola J P B 2010 Lipid Droplets in Inflammation and Cancer *Prostag. Leukotr. Ess.* **82** 243–50
- [4] Herber D L *et al.* 2010 Lipid accumulation and dendritic cell dysfunction in cancer *Nat Med* **16** 880–
- [5] Sun T P *et al.* 2010 The use of bioimpedance in the detection/screening of tongue cancer *Cancer Epidemiol.* **34** 207–11
- [6] Kerhet A, Raffetto M, Boni A and Massa A 2006 A SVM approach to microwave breast cancer detection *Eng. Appl. Artif. Intel.* **19** 807–18
- [7] Bellorofonte C, Vedruccio C, Tombolini P, Ruoppolo M and Tubaro A 2005 Non-invasive detection of Prostate Cancer by Electromagnetic Interaction. *Eur. Urol.* **47** 29–37
- [8] Markx G H and Davey C L 1999 The dielectric properties of cell at radiofrequencies: Applications in biotechnology *Enzyme Microb. Tech.* **25** 161–71
- [9] Montgomery G A 1947 Techniques of microwave measurements (New York, NY: McGraw-Hill)
- [10] Murthy V R K and Raman R 1989 A method for the evaluation of microwave dielectric and magnetic parameters using rectangular cavity perturbation technique *Solid State Commun.* **70** 847–50
- [11] Sheen J 2005 Study of microwave dielectric properties measurements by various resonance techniques *Measurement* **37** 123–30
- [12] Fang Y, Ferrie A M, Fontaine N H, Mauro J and Balakrishnan J 2006 Resonant Waveguide Grating Biosensor for Living Cell Sensing *Biophys. J.* **91** 1925–40
- [13] Rudan-Tasič D and Klofutar C 1999 Characteristics of vegetable oils of some Slovene manufactures *Acta Chim.* **46** 511–21